neptune[™]
reclining bath lift

Operating Instructions

 mountway[™]
bringing independence to life

OPERATING INSTRUCTIONS AND PARTS LIST

Please read the following instructions carefully before using your Neptune Reclining Bath Lift

This bath lift is designed and manufactured to assist individuals who experience difficulties getting in and out of the tub. It should not be used for any other purpose. If you have a helper in the bathroom, ensure they also read these instructions before use. If you are unable to get out of the tub unaided, do not bathe without someone on call to help in an emergency.

The Neptune Reclining Bath Lift is powered by a small rechargeable battery fitted into the hand control. The battery must be inserted into the hand control to power the bath lift.

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Parts list / Liste des pièces

15. Accessories / Accessoires 12/25
16. Parts list / Liste des pièces 13/26

1. Unpacking the Neptune Reclining Bath Lift


Once you have opened the packaging, ensure that the Neptune Reclining Bath Lift has not been damaged in any way during transit. If it has, contact your supplier immediately.

In the box you will find:

- A Neptune Reclining Bath Lift, complete with two side flaps and hand control attached. (see fig 1)
- A battery stick (see fig 2)
- A battery charger (plug and cradle) (see fig 2)

NOTE:

- Ensure that the battery stick is placed on charge for at least 10 hours prior to first use of the bath lift.
- Never transport the Neptune Reclining Bath Lift with the battery stick connected to the hand control.

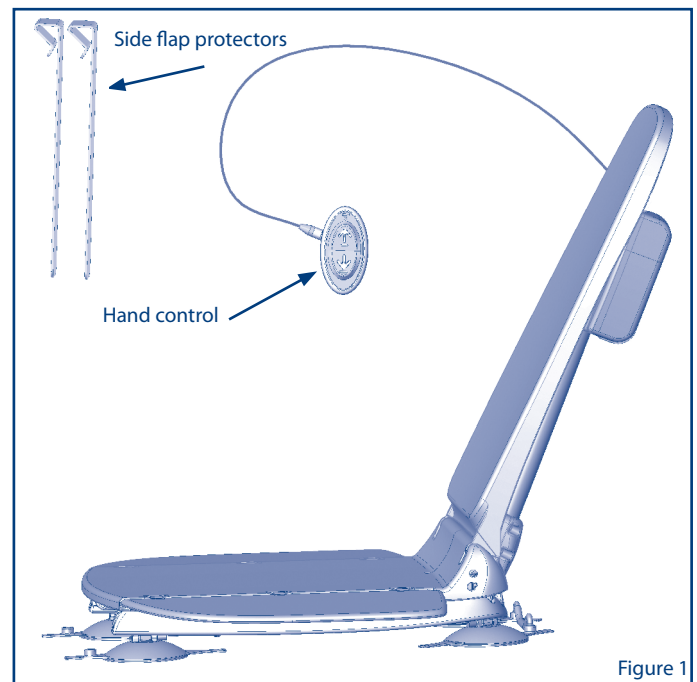2. Recharging the battery


Ensure that the battery stick is placed on charge for at least 10 hours prior to first operation or after long periods of non-use, typically 7-10 days.

- Slide the battery stick into the track on the charger cradle until it clicks into place (see fig 2).
- Plug the charger into the mains electricity supply outside the bathroom and switch on.

With the battery stick inserted into the charger cradle and the electricity turned on, the red charge light appears on the plug. After approximately 10 hours (the time required to charge a flat battery), the red charge light goes off and the green trickle charge light appears. The green trickle charge light indicates that the battery is fully charged and is now ready for use. The battery can also remain continuously on trickle charge until the next bath.

NOTE:

- Do not connect the battery stick to the hand control until the bath lift is fully assembled in the tub.
- When the bath lift is fully assembled in the tub, the battery stick should remain connected to the hand control between each use. It should only be disconnected when the FLASHING GREEN light appears on the hand control to indicate recharging is required.

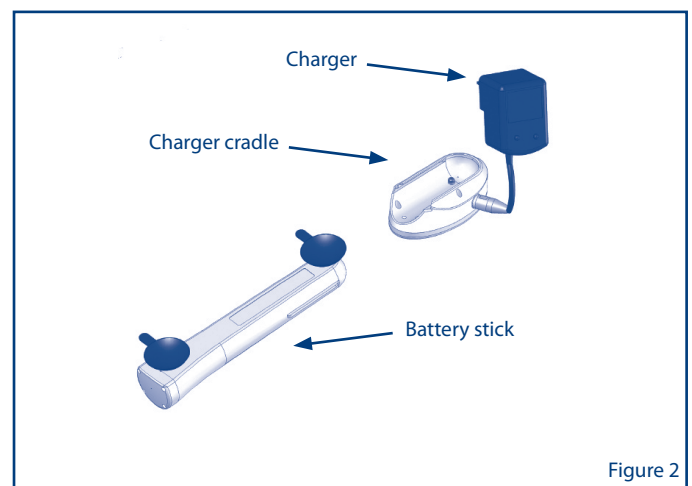3. Preparing the tub

- i) Remove any non-slip mats from the tub.
- ii) Ensure that the bottom of the tub is clean and smooth and has no bath oils or tub cleaner residue that will prevent the bath lift's suction feet from adhering.
- iii) The Neptune Reclining Bath Lift suction feet may not adhere to tubs that have been re-enamelled or have a rippled or textured surface. For some tubs with rippled surfaces, the front suction feet can be relocated to the middle of the base frame away from the textured area. A centre suction foot fixture kit (Part No: AB228) is available on request.
- iv) Ensure no obstructions will impede the lifting movement of the bath lift in the tub e.g. tub handles, soap dishes. The side flap protectors provided should be fitted to tub handles (see fig 3) to prevent the side flaps from catching.

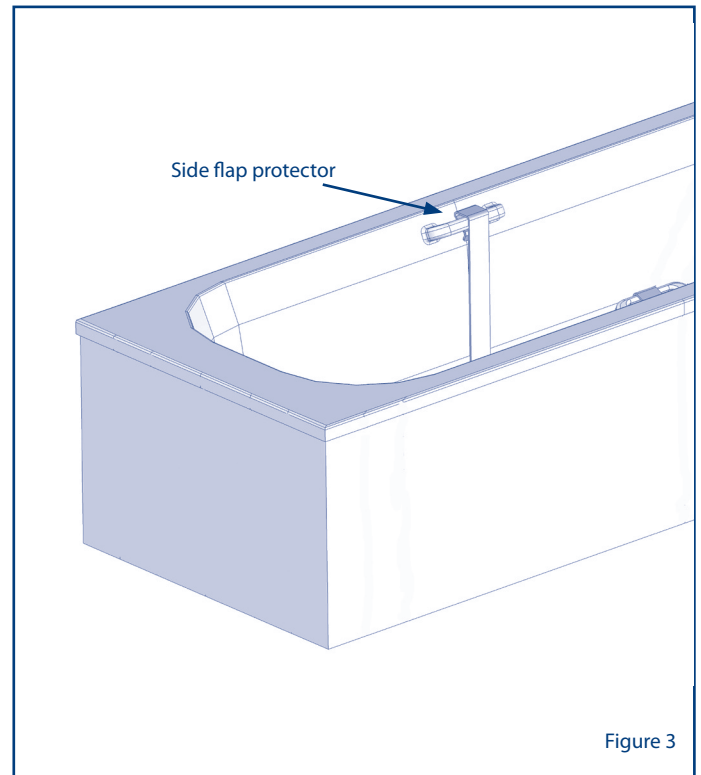
Figure 3

4. Installing the Neptune Reclining Bath Lift

We recommend that the bath lift is always separated and transported in two lighter, individual sections (i.e. backrest assembly and seat unit).

- i) If your Neptune Reclining Bath Lift is supplied with comfort covers (optional), fit them to the bath lift at this stage.
- ii) Place the backrest carefully on the floor when not attached to the seat unit. The actuator forms a useful handgrip for transporting the backrest (see fig 4)
- iii) Pick up the bath lift seat unit. The rear of the base frame forms a useful handhold for this purpose (see fig 5).

NOTE:





- Never use the release handle on the rear suction feet to transport the seat base

The bath lift can be carried as a complete unit in its flat folded position i.e. with the backrest folded onto the seat. (see fig 6)

NOTE:

- The backrest is not locked to the seat when in the folded position and care should be taken when handling as a complete unit.

- iv) Place the seat unit in the bottom of the bath at the opposite end to the taps.
- v) Place as far back as possible while ensuring there is enough clearance at the rear of the tub for the backrest to fully recline. The backrest should not rest on the tub or wall when reclined.
- vi) Press firmly on the seat, ensuring all four suction feet are attached to the bath surface. All suction feet swivel and will adhere to the contours of the tub as required.
- vii) To fit the backrest to the seat unit, place it face down onto the seat, whilst simultaneously aligning the lower two hinges with the slots located at the rear of the seat unit (see Fig 7).
- viii) Hold the top of the backrest and lift upwards until the backrest engages into the reclined position.

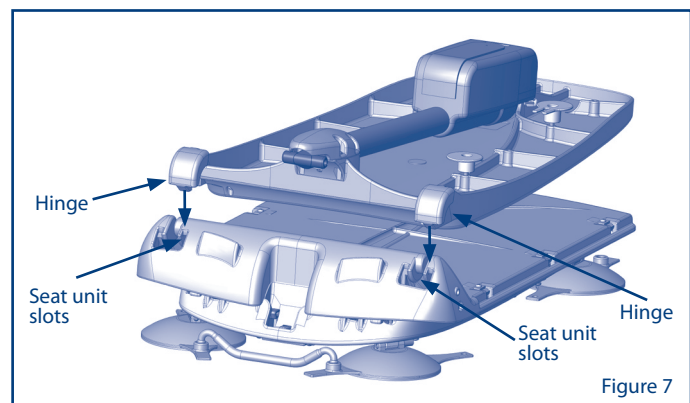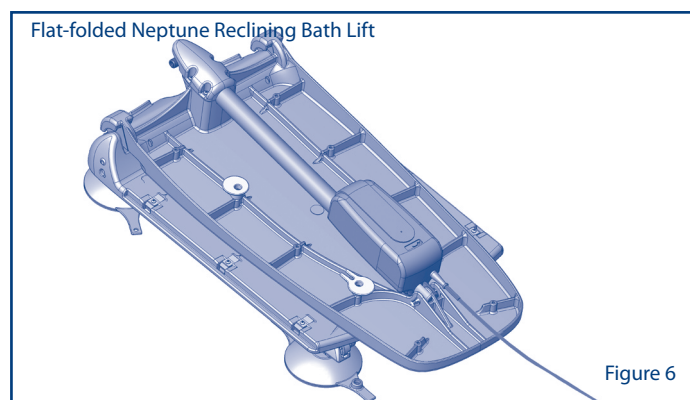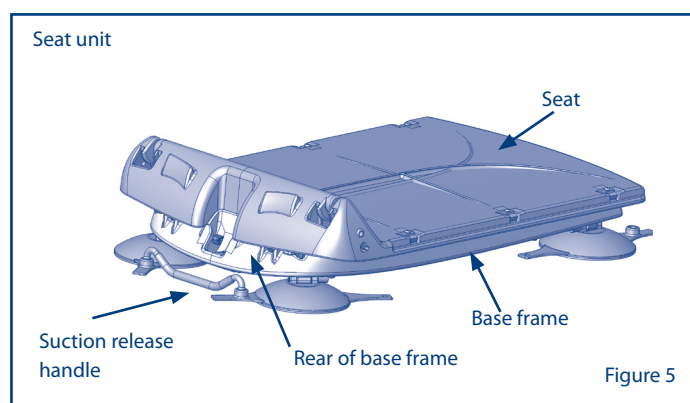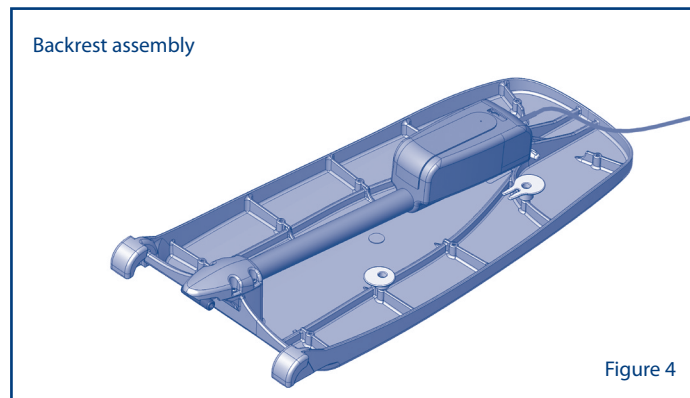NOTE:

- Ensure the actuator end is fully retracted before assembling the backrest to the seat unit.
- Ensure both hinges are fully located in the seat slots.

- ix) Unfold the two side flaps so they rest against the sides of the tub. (The side flaps are fitted to the bath lift for ease of transfer when entering and leaving the tub. They are also there for your safety, as they prevent finger entrapment).
- x) The battery stick can now be connected to the hand control on the bath lift.
- xi) Fill the tub with water; the Neptune Reclining Bath Lift is now ready for use.

5. The hand control

Ensure that the battery stick is placed on charge for at least 10 hours prior to first operation or after long periods of non-use (see section 2 "Recharging the battery").


Slide the battery stick into the track of the hand control until it clicks into place.

The hand control has a battery status indicator, which shows how much charge is in the battery (see fig 9). By observing this indicator you can check the battery status whilst bathing.

- If the light is continuous GREEN, there is sufficient charge for you to have at least one bath.
- If the GREEN light starts to FLASH, you must recharge your battery after this bath (see section 2 "Recharging the battery"). Don't worry if the light starts to flash while you are bathing – there will be enough charge for you to complete your bath.
- If the light turns RED, there will always be enough charge to raise you back to the top of the tub. You must recharge your battery IMMEDIATELY.
- The RED light signals that a safety mechanism has been activated, which prevents the bath lift from lowering you to the bottom of the tub if there is insufficient charge to bring you back up again.

NOTE:


- You should NEVER start a bath when the hand control light is RED.

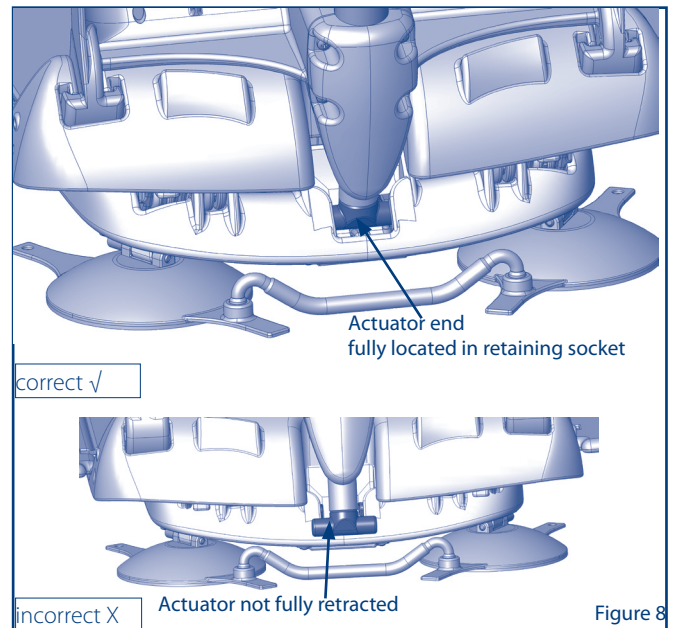- If the battery status indicator was continuous GREEN at the end of your last bath, there will be sufficient charge to continue with your next bath, provided this occurs within seven days. If in doubt, recharge the battery before use.

NOTE:

- The Neptune Reclining Bath Lift is not designed for continuous operation. Allow at least 15 minutes between baths.

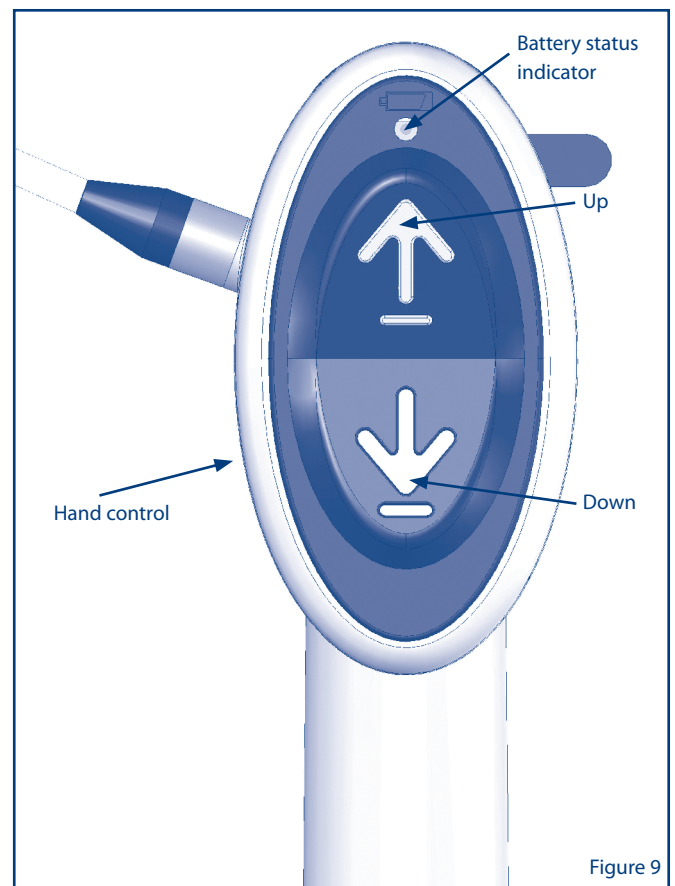6. Using the Neptune Reclining Bath Lift

- i) To raise the bath lift, press the UP button on the hand control (see fig 9) and the bath lift will rise.

The Neptune Reclining Bath Lift is a powerful lifting device and it is essential that the side flaps do not catch under handrails, soap dishes, tile edges, vanity shelves or other obstructions when rising. If the side flaps are caught, immediately lower the bath lift by pressing the DOWN button on the hand control (see fig 9). Pull each side flap free from the obstruction, then resume the lifting operation. (The side flap protectors provided should allow side flaps to glide smoothly over tub handles).

As the Neptune Reclining Bath Lift reaches the height of the tub, the side flap will fold out to bridge the gap between the seat and the side of the tub.

- ii) With the seat and side flap level, sit on the flap and move onto the seat.
- iii) Turn as you lift each leg over the side of the tub then make yourself comfortable on the seat. You may find it easier to transfer if you wet the seat and side flaps first. A slide and swivel disc (Part No: AB238) is available on request, should you still find transferring difficult.
- iv) Extend your feet in front of you and keep your hands above the seat.

NOTE:

- Ensure your legs and hands are clear of the mechanism underneath the seat.

- v) Using the hand control, press the DOWN button (see fig 9) and you will be lowered gently. When the button is released, the bath lift will stop immediately and will remain at that position until the button is pressed again.
- vi) Continue to press the DOWN button until the bath lift reaches its lowest position at the bottom of the tub. By continuing to press the DOWN button, the backrest will begin to recline. Release the DOWN button when you find your desired position. By maintaining pressure on the DOWN button, the backrest will stop at its fully reclined position. To reduce the backrest angle, press the UP button.

- vii) To raise the bath lift, press the UP button. The backrest will move into the upright position and the bath lift will begin to rise.

- viii) Once the side flap is level with the top of the tub, remove your finger from the UP button and move sideways onto the side flap.

- ix) Turn whilst lifting each leg out of the tub.

NOTE:

- If the RED light (indicating low charge) appears on the hand control whilst descending, a safety feature will be activated which will prevent you from going down any further. Press the UP button and the bath lift will raise you safely to the top of the tub.
 - When the Neptune Reclining Bath Lift reaches its maximum limit of travel (up or down), there will be a short delay before it will operate again.
- If the hand control and battery stick fall into the water, there is NO danger as both are waterproof.
- In an emergency, sliding the battery stick out of the hand control can instantly terminate all movement.

7. Removing the Neptune Reclining Bath Lift from the tub

- i) Lower the bath lift until the backrest fully reclines.
- ii) Fold both side flaps onto the seat.

NOTE:

- Always remove the battery stick from the hand control before separating the backrest from the seat unit.

- iii) Drain water from tub, then fold the backrest forward until it lies flat on the seat.
- iv) Wind the hand control cable around the two cable winders on the rear of the backrest, then attach to the retainer clip (see fig 10).
- v) Hold the actuator and lift the backrest upwards and clear of the seat unit. Place carefully on the floor out of the way.
- vi) Pull the handle attached to the rear suction feet and the rear feet will release (See fig 11).
- vii) Reposition your hand to hold the rear of the base

NOTE:

- Never use the suction feet release handle to lift or transport the seat unit.

- viii) Continue to tilt forward and the front suction feet will automatically release.
- ix) Now firmly grip the base frame and lift out of the tub.

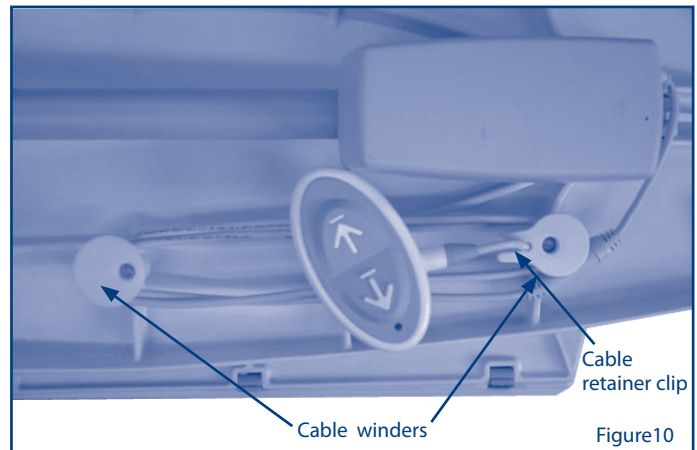
Figure10

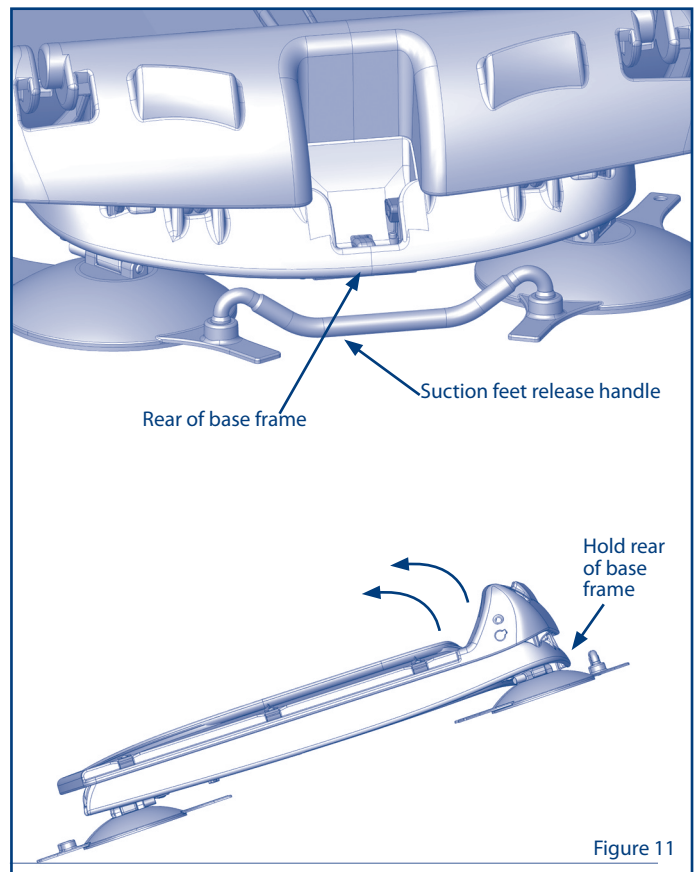
Figure 11

8. Storing the Neptune Reclining Bath Lift

When the bath lift is not installed in the tub, it may be stored with the backrest folded flat onto the seat. We recommend that the hand control cable is wound around the cable winders on the rear of the backrest and then attached to the retainer clip (see fig 10).

NOTE:

- The backrest will not be locked to the seat in this position.
- Do not attach the battery stick to the hand control when the backrest is folded flat.
- Only operate the bath lift when the backrest is assembled and locked to the seat unit.

9. Cleaning

Clean the Neptune Reclining Bath Lift as frequently as you clean your tub, using the same non-abrasive cleaners. Should you wish to disinfect it, use a soft brush or cloth and wipe down with a proprietary sterilising solution or anti-bacterial cleaner. Cleaning agents should be used in accordance with the manufacturer's guidelines. Rinse the product with clean water. Wipe and dry using a clean cloth. Bath oils and other additives are generally easily removed with washing-up liquid. Do not allow oils or other additives to build up.

It is preferable to remove oily deposits frequently as, in the course of time, they may lead to degradation of the suction feet. SOLVENT-based cleaners, such as white spirit, should NOT be used.

The bath lift can be washed in an automated washing machine (to a maximum temperature of 163° F for 3 minutes), however the battery MUST be removed and washed by hand.

If your bath lift is fitted with comfort covers, ensure you remove them on a regular basis to clean the bath lift surfaces.

10. Maintenance and servicing

The Neptune Reclining Bath Lift is designed to be maintenance-free for many years.

However, your Mountway authorised dealer should inspect it annually. The dealer will inspect the bath lift to ensure that there are no loose fittings or other evidence of mechanical damage.

There are no user serviceable parts in the Neptune Reclining Bath Lift. Repairs and servicing must be entrusted only to your Neptune Reclining Bath Lift dealer.

11. Warranty

This product is sold with the benefit that, in the event of any defect of manufacture or material appearing on the frame within 60 months of the date of receipt of the product, the defect will be rectified free of charge, provided that:

- i) Reasonable evidence is provided (e.g. purchase invoice or delivery note) that the product was purchased not more than 60 months prior to the date of claim; and
- ii) The defect is not attributable to accidental damage (either in transit or otherwise), fair wear and tear, wilful damage, misuse or unauthorised repair.

The battery is guaranteed for 24 months to power (when fully charged) one bathing cycle at full load (140kg). Other consumables such as suction feet and charger are also guaranteed for 24 months.

Consumers' statutory rights are not affected.

The bathing cycle:

For one bath, it is assumed that the Neptune Reclining Bath lift will be raised from the bottom to the height of the tub and then lowered into the water supporting a person of up to 140kg. After bathing, the bath lift will raise the person to the height of the tub then be lowered to the bottom of the tub for storage.

12. Precautions

- Do not use the Neptune Reclining Bath Lift for uses other than to assist with getting in and out of the tub and do not use it in combination with other lifting or moving devices.
- To ensure that the hand control is not accidentally damaged, do not allow the cable to trail when transporting or storing. Wind the cable around the two cable winders on the rear of the backrest and then attach to the retainer clip.
- Always check that the bath lift's suction feet are all firmly attached to the tub before use.
- When assembling the backrest, ensure that both hinges are located in the slots on the seat unit and the actuator end is fully retracted and locates correctly into the base frame socket (see fig 8).
- As the Neptune Reclining Bath Lift rises in the tub, ensure that the side flaps do not become caught under any obstructions, such as tub handles, soap dishes or tile edges. Ensure tubs with integral handles are fitted with the side flap protectors provided.
- Before transferring onto the bath lift, ensure that the side flaps are level with and supported by the side of the tub.
- Ensure that no part of the body can become trapped when transferring onto the seat or when raising and lowering yourself on the Neptune Reclining Bath Lift.
- The battery stick is for use solely with Mountway bath lifts and cannot be used to operate any other product.
- The battery stick should not be connected to the hand control and operated without the backrest being fully assembled and locked to the seat unit.
- Ensure that no metal objects touch the contacts on the battery stick, as this will damage the battery.
- The bath lift should not be used under a constant flow of water i.e. running tap. The actuator must not be totally submerged in water.
- The hand control and battery stick are waterproof. If the hand control or battery stick has been immersed in water, they remain perfectly safe but should be wiped dry.
- The battery stick must not be subjected to high temperature cleaning.
- Use only the Mountway charger to charge the Neptune Reclining Bath Lift battery and for no other purpose. Never use any other charger. Never charge in the bathroom or any humid area. Always charge in a dry area, ensuring your hands are dry first.
- When the bath lift is used communally, we recommend cleaning the seat with a suitable sterilising solution or anti-bacterial cleaner after each use. Rinse the seat thoroughly with water before use. Any comfort covers fitted should also be regularly removed and washed.
- DISPOSAL: The plastic, metal parts, packaging materials and suction feet are recyclable but the battery stick, which contains the Nickel Metal Hydride (NiMH) battery, should be disposed of only by returning it to your supplier, the manufacturer or an authorised disposal point.

13. Symbols

To comply with international standards for medical equipment, the following symbols are found on the Neptune Reclining Bath Lift labels:

	This means that only direct current from a small battery is used to power your Neptune. Reclining Bath Lift There is NO connection to mains electricity when in use, so there is no risk of electric shock.
	This means that the Neptune Reclining Bath Lift is categorised as type B medical equipment and meets international standards for electrical safety.
IP66	This means that the bath lift is sealed to the international standard against the effects of water spray.
IP67	This means that the battery stick is sealed to the international standard against the effects of immersion in water to a depth of 1 metre.
IP68	This means that the hand control is protected against long periods of immersion in water under pressure.
	This symbol is found on the hand control and means the light represents the battery status.
	This means PLEASE READ THE INSTRUCTIONS BEFORE INSTALLATION OR USE.
	Do not bin/dispose.


14. Technical data

The Neptune Reclining Bath Lift seat is made of glass-filled polypropylene / nylon. The hand control and battery are made from ABS. The hand control button cover is made of Thermoplastic Elastomer.

The bath lift is in the official category of “internally powered equipment” and is not suitable for use in the presence of flammable anaesthetic mixture with air or nitrous oxide.

The equipment conforms to the requirements of:

- The European Council Directive 89/336/EEC relating to Electromagnetic Compatibility.
- The European Council Directive 93/42/EEC relating to Medical Devices.
- The European Council Directive 73/23/EEC relating to Electrical Safety.

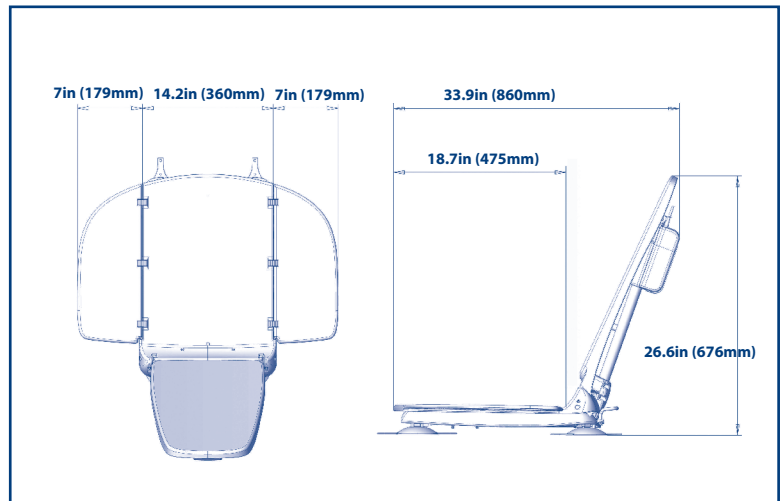Neptune Reclining Bath Lift

Lifting capacity (person)	308 lbs (140 kg)
Net weight (excluding battery)	21.2 lbs (9.6 kg)
Seat base weight	13.2 lbs (6kg)
Backrest weight	7.9 lbs (3.6kg)
Seat height UP	17.3 in (440 mm)
Seat height DOWN	2.9 in (75mm)
Backrest height	26.77 in (680mm)
Power input	14.4V dc nominal 7A max
Duty cycle	One bath per 15 minutes
Degree of protection	IP66
Angle of recline	35°
Hand control (battery)	
Input	17.4V dc nominal 200mA
Output	14.4V dc nominal 7A max
Degree of protection	IP 68
Battery stick	IP67
Battery charger	
Input	6VA 110V ac
Output	18V dc nominal 200mA

Operating and charging

Temperatures	
Max	104°F (40°C)
Min	0°C (32°F)

Shipping and storage

Weight (packaged)	25.4 lbs (11.5 kg)
Shipping volume	2.48ft³ (0.07 m³)
Temperatures	
Max	104°F (40 °C)
Min	14°F (-10 °C)
Humidity	
Max	90% RH
Min	10% RH
Operating Humidity	30 - 75%
Pressure	
Max	125 kPa (18.4 psi)
Min	73 k Pa (10.7 psi)

This product is suitable for transport in pressurised aircraft holds.

15 Accessories

Neptune Reclining Bath Lift can be supplied with a wide range of optional accessories:

Item

1. Pommel (Part No. AB248)

– Separates legs for maintaining correct posture or improving hygiene access.

2. Knee wedge (Part No. AB249)

– Supports legs, maintains correct posture or prevents sliding.

3. Soft cover set (back & seat) (Part No. AB224)

– Increases comfort.

4. Cushioned cover set (back & seat) (Part No. AB225)

– Provides even greater comfort.

5. Cushioned cover set paediatric (back & seat) (Part No. AB246)

– For use with pommel and knee wedge.

6. Thoracic harness adult inc fixture kit (Part No. AB232)

– Helps maintain upper body posture or prevents slipping.

7. Thoracic harness child inc fixture kit (Part No. AB231)

– Helps maintain upper body posture or prevents slipping.

8. Pelvic harness adult (Part No. AB230)

– Helps maintain body posture or prevents slipping.

9. Pelvic harness child (Part No. AB229)

– Helps maintain body posture or prevents slipping.

10. Slide and swivel aid (Part No. AB247)

– Facilitates transfer on and off the seat.

11. Height adapter set (25mm) (Part No. AB213)

– Fits to the suction feet to increase height of bath lift in exceptionally deep baths.


12. Center suction cup fixture kit (Part No. AB228)

– Attaches front suckers to middle of base frame in baths with rippled surfaces.

13. Carry bag (Part No. AB250)

– Protects the bath lift while travelling and facilitates transport.

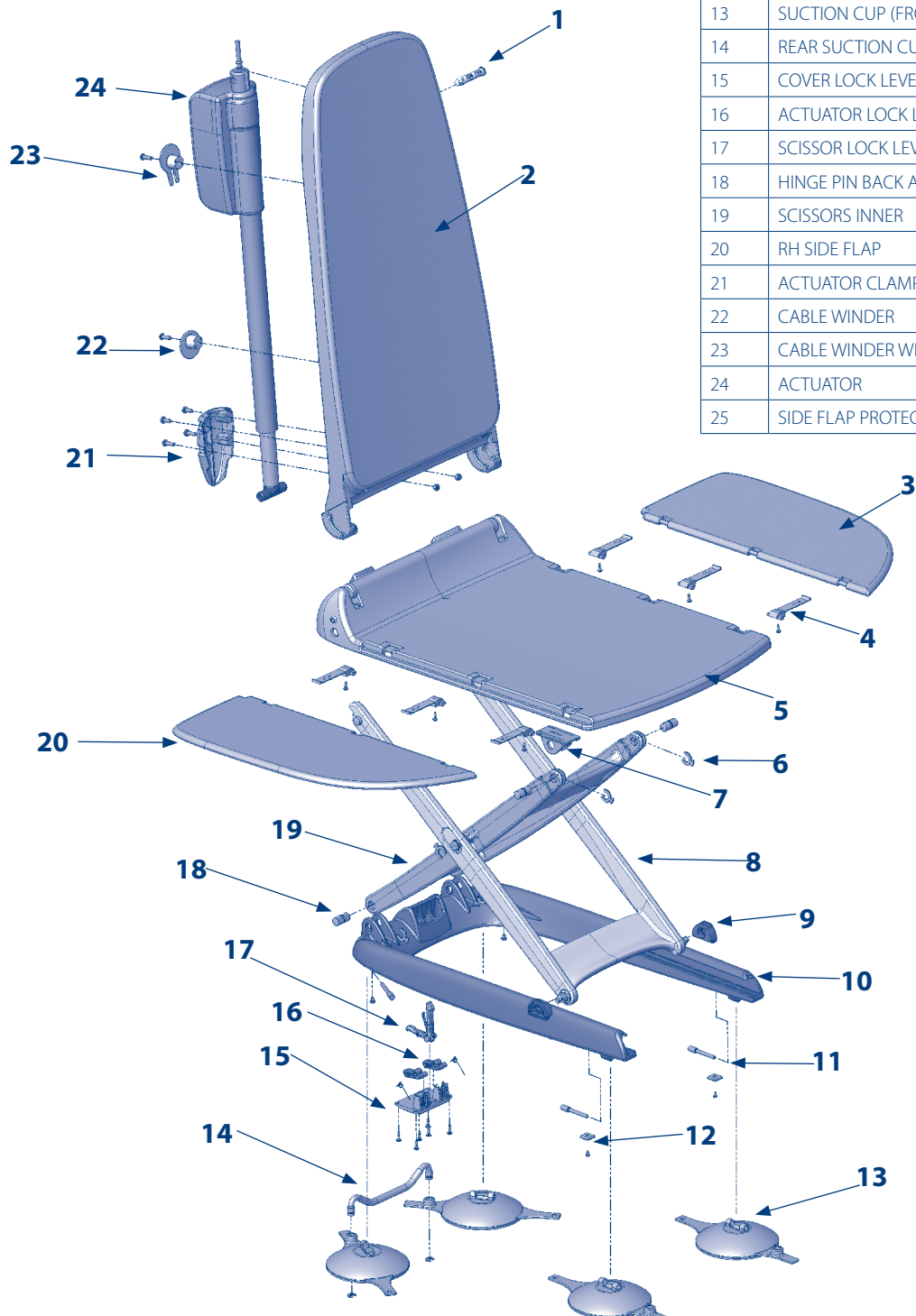ITEM NO.	DESCRIPTION
1	ACTUATOR SNAP PIN
2	BACK REST
3	LH SIDE FLAP
4	HINGE
5	SEAT
6	SNAP WASHER
7	SLIDE BLOCK UPPER
8	SCISSORS OUTER
9	SLIDE BLOCK LOWER
10	BASE FRAME
11	SUCKER SWIVEL PIN- FRONT AND REAR
12	RECTANGULAR WASHER
13	SUCTION CUP (FRONT)
14	REAR SUCTION CUP ASSEMBLY
15	COVER LOCK LEVERS
16	ACTUATOR LOCK LEVER
17	SCISSOR LOCK LEVER
18	HINGE PIN BACK AND FRONT
19	SCISSORS INNER
20	RH SIDE FLAP
21	ACTUATOR CLAMP
22	CABLE WINDER
23	CABLE WINDER WITH RETAINER CLIP
24	ACTUATOR
25	SIDE FLAP PROTECTOR (NOT SHOWN)

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Neptune is manufactured by Mountway Limited